Sunbeam

Model 884

Extended Back Massaging Renu™ Tension Relieving Heat Therapy Pad

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PRODUCT IS NOT TO BE USED ON OR BY AN INVALID, A SLEEPING OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, OR A PERSON WITH POOR BLOOD CIRCULATION.
5. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
6. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PRODUCT FREQUENTLY.
7. DO NOT USE IN AN OXYGEN ATMOSPHERE.
8. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PRODUCT IN PLACE.
9. DO NOT SIT ON OR AGAINST OR CRUSH PRODUCT—AVOID SHARP FOLDS.
10. NEVER PULL THIS PRODUCT BY THE SUPPLY CORD.
11. DO NOT USE THE CORD AS A HANDLE.
12. THIS UNIT SHOULD NOT BE USED BY CHILDREN WITHOUT ADULT SUPERVISION.
13. PLACE PRODUCT ON TOP OF AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.
14. USE THIS PRODUCT ONLY ON A 120 VOLT AC CIRCUIT. UNPLUG WHEN NOT IN USE.
15. DO NOT TAMPER WITH THIS PRODUCT IN ANY WAY.
16. THERE ARE NO USER-SERVICEABLE PARTS. IF FOR ANY REASON THIS PRODUCT DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CONSUMER SERVICE INFORMATION.
17. DO NOT USE THIS PRODUCT WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS. SKIN BURNS COULD RESULT.
18. DO NOT BEND OR PINCH CORD AND DO NOT USE IF CORD IS DAMAGED.
19. NEVER LEAVE THE PRODUCT UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
20. LOOP CORD LOOSELY WHEN STORING. TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
21. NEVER USE PAD WITHOUT COVER IN PLACE.
22. CAREFULLY EXAMINE INNER COVER BEFORE EACH USE. DISCARD THE PAD IF INNER COVERING SHOWS ANY SIGN OF DETERIORATION.
23. SAVE THESE INSTRUCTIONS.

THIS APPLIANCE HAS A POLARIZED PLUG (ONE BLADE IS WIDER THAN THE OTHER). TO REDUCE THE RISK OF ELECTRIC SHOCK, THIS PLUG IS INTENDED TO FIT IN A POLARIZED OUTLET ONLY ONE WAY. IF THE PLUG DOES NOT FIT FULLY IN THE OUTLET, REVERSE THE PLUG. IF IT STILL DOES NOT FIT, CONTACT A QUALIFIED ELECTRICIAN TO INSTALL THE PROPER OUTLET.
DO NOT ATTEMPT TO MODIFY THIS PLUG OR DEFEAT THIS SAFETY FEATURE IN ANY WAY.

Extended Back Massaging Renu™ Tension Relieving Heat Therapy Pad is specifically designed to fit the neck, shoulders & upper back area, where users experience the most physical signs of tension and stress. Drape the Extended Back Massaging Renu™ Heat Therapy Pad comfortably on your shoulders making sure that the neck flap is properly in place. Weighted edges and adjustable snap enclosures help to secure the Extended Back Massaging Renu™ Heat Therapy Pad in place.

HOW TO USE YOUR EXTENDED BACK MASSAGING RENU™ TENSION RELIEVING HEAT THERAPY PAD

The Extended Back Massaging Renu™ Tension Relieving Heat Therapy Pad has 2 heat settings and 4 massage settings which can be used simultaneously or individually for customized therapy (Figure 1).

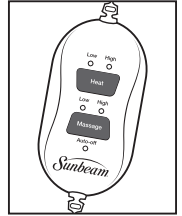
Figure 1

1. Plug into 120V AC electrical outlet.
2. To turn on select either Low/ High Heat for heat therapy, or Low/ High Massage for massage therapy.
NOTE: The 4 Massage settings consists of the following speeds which can be used along with the heat therapy options: (1) Low vibration - consistent, (2) Low vibration - pulsing, (3) High vibration - consistent, (4) High vibration - pulsing.
3. Select desired heat and/ or massage level. **To turn OFF, you must scroll through the remaining heat and/ or massage settings by pressing the buttons sequentially.**
4. Allow the Extended Back Massaging Renu™ Heat Therapy Pad to preheat for several minutes. The Extended Back Massaging Renu™ Heat Therapy Pad has been designed to meet the maximum heat output according to UL Standards.

NOTE: A red blinking light indicates an error has occurred. Unplug from outlet, verify physical condition of the cord and function of AC outlet. If blinking light persists, discontinue use and contact Consumer Services.

AUTO-OFF FEATURE

Your Extended Back Massaging Renu™ Heat Therapy Pad controller includes an Auto-Off feature. This feature is designed to conserve energy by automatically turning off after 2 hours of continuous use when using the Heat with Massage option, and a 30-minute Auto-Off feature when using only the Massage option. The red Auto-Off light will blink, indicating that this feature has activated. **NOTE:** After a power outage, the Auto-Off feature automatically resets the timer.

WASHING INSTRUCTIONS

**DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER!
DO NOT WET INNER PAD CONNECTED WITH CORD! ONLY WASH REMOVABLE OUTER COVER!**

1. **CAUTION:** Disconnect the Extended Back Massaging Renu™ Heat Therapy Pad from outlet.
2. Remove the outer cover from the Extended Back Massaging Renu™ Heat Therapy Pad.
3. Machine wash the cover ONLY in cold water on gentle cycle, or hand wash.
4. Tumble dry on low heat or line dry.

2-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of two years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or in any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition.

Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise.

JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A. - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-435-1250 and a convenient service center address will be provided to you.

In Canada - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-667-8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Herford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

©2016 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved. Distributed by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions, Boca Raton, Florida 33431.
Printed in China. DW-042116

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