

TREATMENT AREAS

The Trinity ELE is designed to stimulate smaller surface areas of the face, such as the lip and eye areas. The instructions below are designed to provide users with the most effective results when followed correctly.

Always ensure that you are correctly holding the device in either a vertical or horizontal position, as required per the specified treatment area. Reference the Treatment Key below as a guide on how to properly position the spheres during treatment.

TREATMENT KEY

Arrows represent both spheres in a gliding action:



Horizontal gliding action.



Vertical gliding action.

Dots represent both spheres in a holding (stationary) action:



Horizontal holding action.



Vertical holding action.

LIPS


1. LOWER LIP: Place both spheres at Position 1 horizontally under the lip, slightly off center, and hold for 1 beep (5 seconds). Move to Position 2 and hold for 1 beep. Complete sequence 3 times.

2. UPPER LIP: Place both spheres at Position 1 horizontally above the upper lip and hold for 1 beep (5 seconds). Move to Position 2 and hold for 1 beep. Complete sequence 3 times.

3. NASOLABIAL FOLD: Place both spheres at Position 1 vertically starting at the upper corner of the mouth and hold for 1 beep (5 seconds). Move to Position 2 and hold for 1 beep. **DO NOT PLACE THE SPHERES DIRECTLY ONTO OR INTO THE NOSE.** Complete sequence 3 times.

EYES

4. LOWER EYE: Place both spheres at Position 1 horizontally at the outer corner of the orbital rim of the eye and glide inward toward inner corner of the eye until the beep sounds (5 seconds). **DO NOT PLACE THE SPHERES DIRECTLY INTO THE TEAR DUCT.** Repeat glide at Position 2 and then Position 3. Complete sequence 3 times.

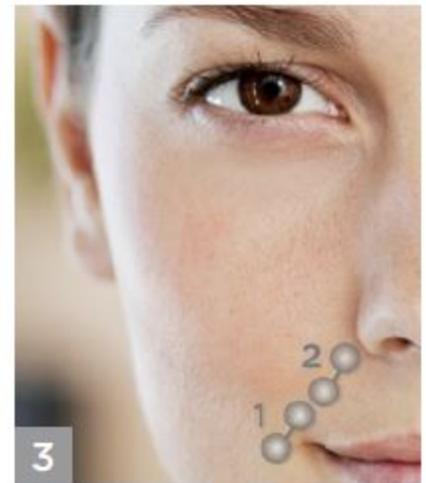


5. UPPER EYE: Place both spheres at Position 1 at the outer corner of the upper orbital rim of the eye and glide inward toward the inner corner of the eye until the beep sounds (5 seconds). **DO NOT PLACE THE SPHERES DIRECTLY ON THE EYELIDS.** Repeat glide at Position 2 and then Position 3. Complete sequence 3 times.



6. EYEBROW: Place both spheres at Position 1 horizontally just below the eyebrow. Lift upward slightly and hold for 1 beep (5 seconds). Move to Position 2 and hold for 1 beep, then to Position 3 and hold for 1 beep. Complete sequence 3 times.



7. IN-BETWEEN BROWS: Place both spheres at Position 1 horizontally in the middle of the eyebrows, slightly off-center. Glide the spheres in an upward motion ending at the center of the forehead until the beep sounds (5 seconds). Complete the glide 3 times.

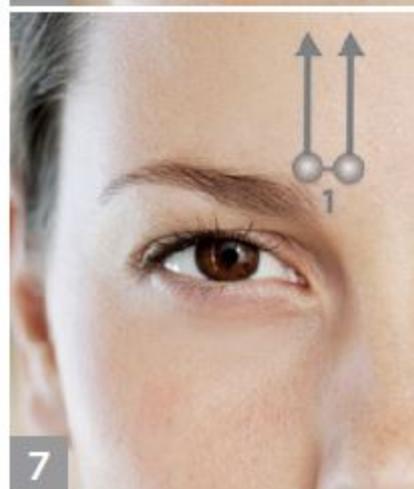

Repeat treatment steps on opposite side of the face.

Your device will automatically shut off after 21 minutes. Keep your NuFACE Trinity device in its charging cradle when not in use. See the “Cleaning” section for proper cleaning instructions. **Note: Do not immerse device or treatment attachment in water at any time.**

TREATMENT FREQUENCY

- Use on targeted areas of concern (3 minutes for the lips, 5 minutes for the eyes).
- For the first 60 days, use 5 times per week.
- After 60 days, use 2-3 times per week for maintenance.
- Allow 24 hours between treatment sessions.

RESULTS

Results may vary. Most users experience instant facial stimulation benefits while others take longer. Optimal results are achieved with regular use.