



Le Mieux

SKIN PERFECTER

rediscover perfect skin

USER GUIDE




3 Easy Steps to Perfect Skin

STEP 1: Prepare




Dampen face and neck with water.

STEP 2: Deep Cleanse, Exfoliate, Extract

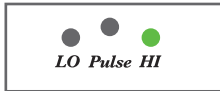


Press **Power** button, then **Mode** button 2 times until HI is green.

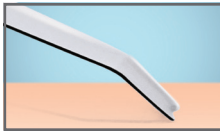


Hold unit with tip of **spatula angled down**. **Keeping face wet**, gently glide spatula outward and upward over face, nose, and neck. Turn off unit. Pat dry face and neck.

STEP 3: Infuse

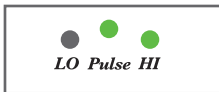


Apply serum or other skincare treatment to face and neck.

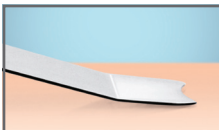

Press **Power** button, then **Mode** button 3 times until HI and Pulse are green. Hold unit with tip of **spatula flat against skin**. Gently glide spatula outward and upward over face, nose, and neck.

Optional: Use circular motions for micro-massage.

Just **3** easy steps, **3** minutes, **3** times a week to rediscover perfect skin.

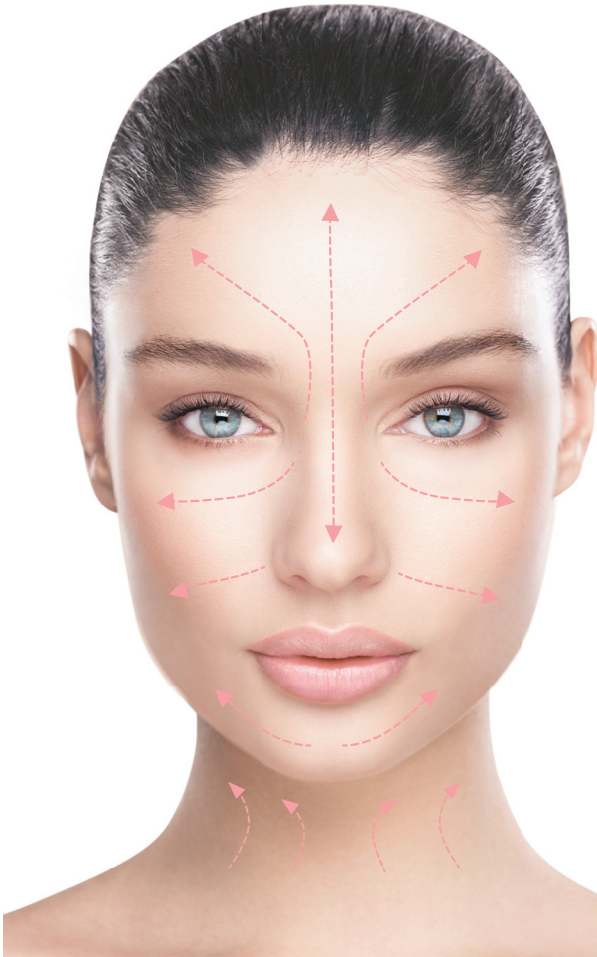

Mode button →



Power button →





4-in-1 Beauty Tool

Beauty really is skin deep.

Deep within is beautiful, younger skin waiting to be revealed.

The Skin Perfecter transcends ordinary cleansing by combining four skin-perfecting treatments into one ultrasonic, handheld beauty tool.

Professional strength ultrasonic waves gently cleanse, exfoliate, extract, and infuse to reveal smoother, brighter, and perfect-looking skin.

Micro-massage action provides 28,000 ultrasonic vibrations per second to gently cleanse, tone, and stimulate the skin while accelerating the absorption of skincare treatments.

Plus, the Skin Perfecter is easy to use and never requires brushes, cartridges, or discs.



4⁻ⁱⁿ⁻¹
beauty
tool

1 CLEANSER

deep cleanses skin and pores of excess oil, makeup residue, and surface impurities.

2 EXFOLIATOR

gently removes cellular build-up using non-abrasive micro-exfoliation.

3 EXTRACTOR

extracts imperfections and reduces the appearance of enlarged pores.

4 INFUSER

accelerates the absorption of serums and other skincare treatments.

Every treatment with the
Skin Perfecter visibly improves:



- Perfect for all skin types
- 4 levels of intensity
- Rechargeable battery
- DC power supply adapter included
- Automatic alarm & shut off after 10 minutes

Parts & Control Panel



Power Supply Adapter



Intensity Indicator	LED Display	Skin Type	Mode
LO	● ● ●	Sensitive	Cleanse
LO + pulse	● ● ●	Sensitive	Infuse
HI	● ● ●	Other	Cleanse
HI + pulse	● ● ●	Other	Infuse

Overview



Treatment Areas:	face, neck
Duration of Treatments:	3 minutes
Frequency of Treatments:	3 times a week

- Wet skin thoroughly when using the Skin Perfecter. If skin becomes dry, rewet with fingers or with a saturated cotton pad.
- Do not rinse the spatula under running water or immerse it in water.
- During treatment, do not press the spatula on the same area of skin for a prolonged period of time.

Directions


I. Prepare

- 1 Charge the Skin Protector, using the enclosed power supply adapter. The Battery Indicator will turn from red to off when the unit is fully charged.
- 2 Wash face with daily cleanser. Pat dry.
- 3 Remove the Skin Perfecter's protective cap. Clean the spatula with water or a cotton pad saturated with alcohol. Dampen face.
Skin must remain wet when using the Skin Perfecter.


II. Cleanse

To deep cleanse, exfoliate, and extract skin imperfections. Switch the unit to the appropriate position, according to skin type.


Sensitive	LO	● ● ●	Press Power button.
Other	HI	● ● ●	Press Power button, then press Mode button 2 times.

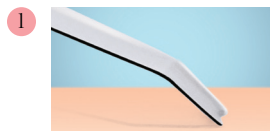


Hold the unit with the tip of the **spatula angled down**.

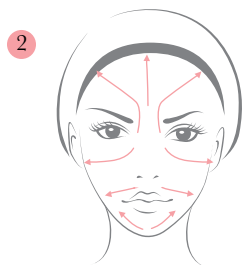

FACE: Keeping the face wet, gently glide the spatula from the center of the face, moving outward.

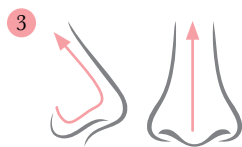

NOSE: Pull the nose to the side. Gently glide the spatula from the edge of the nose toward the center. Repeat process on the other side of the nose. Finish by placing the curved end of the spatula at the bottom of the nose, and glide the spatula upward.

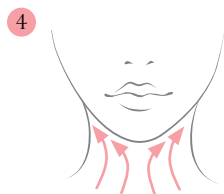

NECK: Dampen neck. Glide the spatula from the bottom of the neck upward, moving toward the bottom of the chin. Turn off the Skin Perfecter and gently pat dry face and neck.



Clean the spatula with water or a cotton pad saturated with alcohol.

III. Infuse

To infuse, tone, and stimulate the skin, first apply skincare product to face and neck. Switch the unit to the appropriate position, according to skin type.


Sensitive	LO + pulse	● ● ●	Press Power button, then press Mode button 1 time.
Other	HI + pulse	● ● ●	Press Power button, then press Mode button 3 times.

1




Hold the unit with the tip of the **spatula flat against the skin**.

2

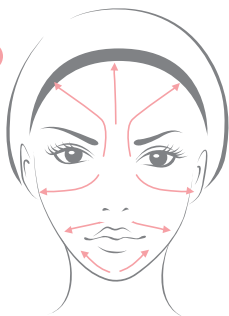


FACE: Gently glide the spatula from the center of the face, moving outward.

3

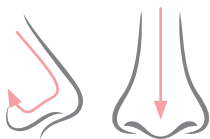

NOSE: Gently glide the spatula from the center of the nose toward the edge. Repeat process on the other side of the nose. Finish by placing the flat end of the spatula at the bridge of the nose, and glide the spatula downward.

4

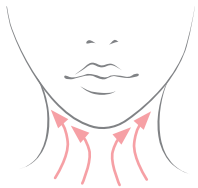

NECK: Glide the spatula from the bottom of the neck upward, moving toward the bottom of the chin. Turn off the Skin Perfecter.

5



Clean the spatula with water or a cotton pad saturated with alcohol.

IV. Finish

Cover with the protective cap. Store in a cool, dry place.

Care & Maintenance

After each use, clean the stainless steel spatula with water or a cotton pad saturated with rubbing alcohol. Cover with the protective cap. Store in a cool, dry place.

Important Safeguards



- Keep out of the reach of children.
- Operate this unit only according to its intended use as described in this manual.
- Do not use this unit if you have serious allergies, infected skin, or a heart condition.
- Do not immerse the unit in water. If the product appears damaged or has fallen into water, discontinue use and return to original place of purchase.
- Do not operate this unit in a room with temperatures higher than 40° C / 104° F or lower than 0° C / 32° F.
- Keep the unit free of lint and hair.

Troubleshooting

After checking the following, if any problem persists, please discontinue use and return to original place of purchase.

Problem	Possible Cause & Solution
Device does not work when powered on.	Battery may be low or empty. Recharge device.
Ultrasonic vibrations are too low, or there is no vibration at all.	Battery may be low or empty. Recharge device. <i>OR</i> The spatula may be deformed or detaching from the handle. Discontinue use and return to original place of purchase.



Specifications

Product Code: LM1401SP

Power Supply Adapter: Input 100-240V AC 50/60Hz; Output 5V DC 2A

Power Consumption: 10 W

Frequency: approx. 28 kHz

Timer: 10 minutes

Battery: Li-Ion 3.7V 900mAh

Battery Life Cycle: approx. 50 hours (300 treatments)

Recharge Time: approx. 3 hours

Battery Usage: approx. 2 hours

Battery Disposal: This device contains lithium batteries. Please follow the regulations where you live for recycling or disposing of electrical devices.

Material of Spatula: SUS

Device: Made in Taiwan

Dimensions: 6.8" x 2" x 0.74" (174 mm x 52 mm x 19 mm)

Weight: 3.88 oz (0.11 kg) without adapter; 7.05 oz (0.2 kg) with adapter



Le Mieux

SkinPerfecter.com

© 2014 Le Mieux Cosmetics, Inc. All rights reserved.

Le Mieux and Le Mieux logo are registered trademarks of Le Mieux Cosmetics, Inc.