

Usage of the LumaRx™ Full-Body Hair Removal System (Indications)

The LumaRx Full-Body Hair Removal System is an over-the-counter device indicated for the removal of unwanted body and/or facial hair in adults with Fitzpatrick skin types I – IV.

The LumaRx Full-Body Hair Removal System is also intended for permanent reduction in unwanted hair.

Permanent hair reduction is defined as the long-term stable reduction in the number of hairs regrowing when measured at 6, 9, and 12 months after the completion of a treatment regimen.

Skin Tone

The LumaRx Full-Body Hair Removal System is intended for use only on light to medium skin tones. To determine if your Skin Tone is suitable to use with LumaRx Full-Body Hair Removal System, compare your skin in the intended treatment areas with the Skin Tone Chart. If it matches the Light to Medium Skin Tones, you can use


LumaRx Full-Body Hair Removal System. If your Skin Tone in the intended treatment area matches the Dark Skin Tones, DO NOT use the LumaRx Full-Body Hair Removal System.


The LumaRx Full-Body Hair Removal System also has a built-in SKIN TONE TESTER that is designed to help prevent flashing on skin that is too dark for treatment. See page 16 for details. Test each different area you wish to treat with the Skin Tone Chart and SKIN TONE TESTER before treating. Some areas are darker than others.

If you have tanned skin in the intended treatment area, DO NOT use the LumaRx Full-Body Hair Removal System in those areas as your skin may be too dark for treatment. You should wait until your tan has completely faded from any intended treatment area before using the LumaRx Full-Body Hair Removal System.

Hair Color

The light energy targets dark coloring in hair, and therefore works best on **DARKER HAIR** (black, dark brown, brown or dark blond). **DO NOT** use the LumaRx Full-Body Hair Removal System on **LIGHTER HAIR** (light blond, red, grey or white hairs or “peach fuzz”). The LumaRx Full-Body Hair Removal System will not work on **LIGHTER HAIR**.

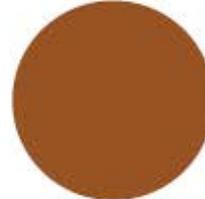SKIN TONE CHART					
					
YES	YES	YES	YES	NO	NO

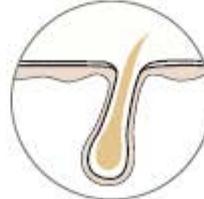HAIR COLOR CHART					
Black	Brown	Dark Blond	Light Blond	Red	Peach Fuzz/ White/Gray
					
YES	YES	YES	NO	NO	NO

What Is the LumaRx™ Full-Body Hair Removal System used for? (Indications)

Treatment Areas


Male and Female Body Sites

The LumaRx Full-Body Hair Removal System is intended for use on male and female body hair on **underarms, bikini line, legs, arms, chest** and **stomach**.

If you have unwanted DARKER HAIR in any of these areas, the illustration below lists the typical number of flashes needed to cover a 1.5" x 3" treatment area. Treat in a GRID-LIKE PATTERN – slight overlaps are OK, but DO NOT flash the same EXACT spot multiple times during one treatment session. Record treatments in the **TREATMENT CARDS** on page 54.

Note: DO NOT physically mark the treatment area with grid pattern as this could cause injury to the skin.

SIZE OF
FLASH
WINDOW

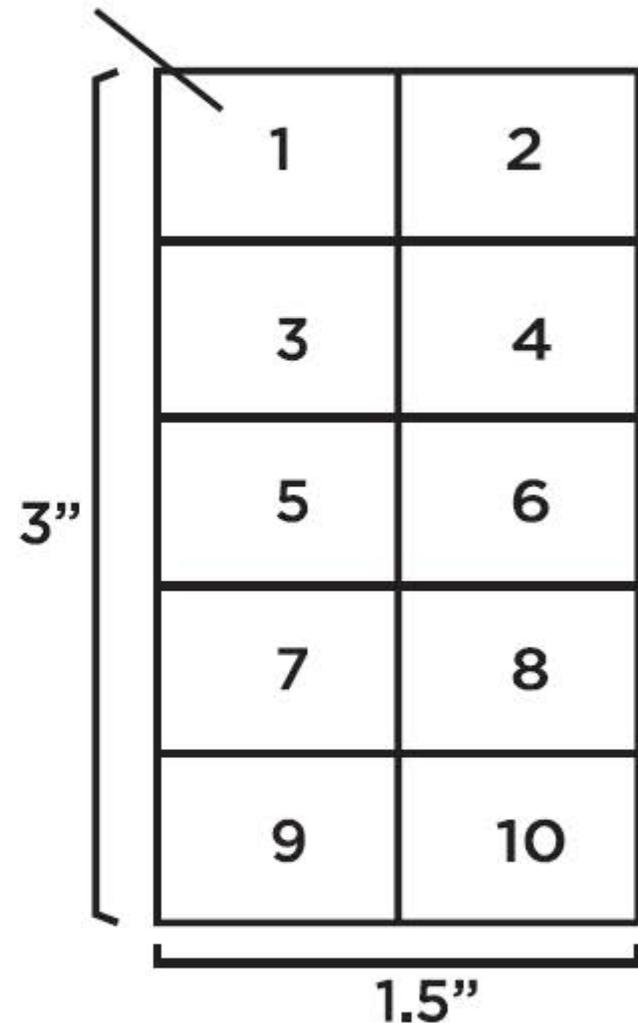


Female Facial Sites

Using the included Facial Cap, the LumaRx Full-Body Hair Removal System is intended for use on female facial hair on the **cheeks, sideburns, chin, neck areas,** and **above the lips** (see illustration). If you have unwanted DARKER HAIR in any of these areas, the illustration below lists the typical number of flashes to treat the entire site (NOTE: for spot treatments of unwanted hair, fewer flashes may be needed.). Record treatments in the TREATMENT CARDS on page 54.

DO NOT use the LumaRx Full-Body Hair Removal System on or around **EYES, EYE LASHES** or **EYE BROWS**. Doing so can cause serious and permanent eye injury.

- Each Sideburn - 5 flashes
- Each Cheek - 15 flashes
- Above lips - 5 flashes
- Chin - 10 flashes
- Neck - 15 flashes

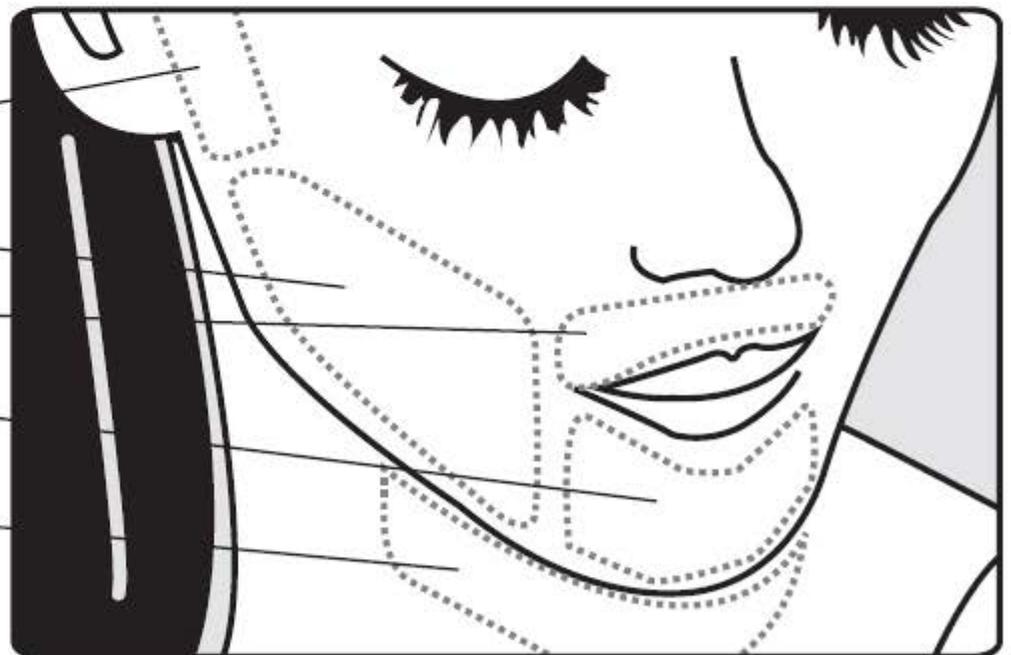

Who shouldn't use the LumaRx™ Full-Body Hair Removal System (Contraindications)

Important Safety Information – Read Before Use!

The LumaRx Full-Body Hair Removal System is not designed for everyone. Please read and consider the information in the following section before use.

- **DO NOT USE** the LumaRx Full-Body Hair Removal System if you have **DARK SKIN**. Darker skin may absorb too much light energy which could cause skin swelling, temporary blisters, or could change the color of your skin. Refer to the Skin Tone Chart on Page 13 to determine if your Skin Tone is suitable for use with the LumaRx Full-Body Hair Removal System. You should test each different area of skin you want to treat. Some areas are darker than others.

- **DO NOT USE** the LumaRx Full-Body Hair Removal System on **TATTOOS, DARK BROWN or BLACK SPOTS** (such as dark freckles, birthmarks, moles, or warts), **LIPS, NIPPLES, GENITALS or AROUND THE ANUS**. These areas have darker skin and may absorb too much light energy which could cause skin swelling, temporary blisters, or could change the color of your skin.
- **DO NOT USE** the LumaRx Full-Body Hair Removal System **ON OR AROUND EYES, EYEBROWS or EYELASHES**. Doing so can cause serious and permanent eye injury.
- **DO NOT USE** the LumaRx Full-Body Hair Removal System anywhere on the **MALE FACE, HEAD OR NECK**.

Warnings

⚠ The LumaRx Full-Body Hair Removal System is a powerful device. As such it should be used with special attention to safety. Please read all warnings and safety precautions before use and follow them strictly when using the LumaRx Full-Body Hair Removal System.

⚠ **DO NOT** flash in an environment where the emission of light could cause an explosion, such as in the presence of gas.

⚠ **KEEP OUT OF REACH OF CHILDREN!** **DO NOT** leave the LumaRx Full-Body Hair Removal System unattended around children as the cords represent a choking hazard. Children should be supervised to ensure that they do not play with the appliance.

⚠ **The LumaRx Full-Body Hair Removal System SHOULD NOT** be used by anyone under the age of 18.

⚠ **DO NOT** use the LumaRx Full-Body Hair Removal System if you have a known **sensitivity to light (photosensitivity)** or are taking medications that make your skin more sensitive to light, as the LumaRx Full-Body Hair Removal System may injure your skin or make existing conditions worse.

Precautions

The LumaRx Full-Body Hair Removal System is a powerful device. As such it should be used with special attention to safety. Please read all warnings and safety precautions before use and follow them strictly when using the LumaRx Full-Body Hair Removal System.

- **DO NOT** use the LumaRx Full-Body Hair Removal System if you have **TANNED SKIN** in the intended treatment area, as your skin may be too dark and treatment in tanned areas could temporarily cause skin swelling, blisters, or could change the color of your skin. Wait until your tan has faded in the treatment area before treating.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on skin with **MAKEUP** applied, as treating in those areas could temporarily cause skin swelling, blisters or could change the color of your skin. Make sure to thoroughly clean the intended treatment area before using the LumaRx Full-Body Hair Removal System.
- **DO NOT** flash the LumaRx Full-Body Hair Removal System on the identical treatment spot multiple times during the same session, as this could temporarily cause skin swelling, blisters, or could change the color of your skin. See instructions on page 32 for more information on how to properly administer treatment.
- **DO NOT** treat any area more than once every other week, as this will not speed up results, but could increase likelihood of skin swelling, blistering or changing the color of your skin.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on damaged skin, including open wounds, cold sores, sunburns, or areas where you have had a skin peel or other skin resurfacing procedure within the last 6-8 weeks.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on irritated skin, including rashes or swollen skin. These conditions can make the skin more sensitive to the light treatment. Wait for

the affected area to heal before using the LumaRx Full-Body Hair Removal System.

- **DO NOT** use the LumaRx Full-Body Hair Removal System if your skin is sensitive to light. In order to determine if you can use the LumaRx Full-Body Hair Removal System, perform a test on a small patch of skin and wait 24 hours to determine if there are any adverse reactions before using the LumaRx Full-Body Hair Removal System more broadly. See steps for testing on a patch of skin on page 30. Contact Customer Service (page 50) for more information or assistance.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on any area where you may someday want your hair to grow back. The results may be permanent and irreversible.
- **DO NOT** flash the LumaRx Full-Body Hair Removal System against any material other than skin.
- **DO NOT** use the LumaRx Full-Body Hair Removal System for any purpose other than hair removal.
- **DO NOT** use the LumaRx Full-Body Hair Removal System if you are or may be pregnant or are breast-feeding. The LumaRx Full-Body Hair Removal System has not been tested on pregnant women. Hormonal changes are common during pregnancy or breast-feeding and could increase skin sensitivity and the risk of skin injury.
- **DO NOT** cover any part of the LumaRx Full-Body Hair Removal System ventilation slots during use, as they help to keep the device cool during operation. Covering the ventilation slots could result in the LumaRx Full-Body Hair Removal System overheating.
- **DO NOT** use the LumaRx Full-Body Hair Removal System if the box packaging or the device appears damaged or if any part of the LumaRx Full-Body Hair Removal System body is cracked or broken. Doing so could increase the risk of skin or eye injury from the flash lamp. Contact Customer Service for assistance (page 50).

Precautions

- **DO NOT** allow the LumaRx Full-Body Hair Removal System to become damp or wet. **DO NOT** use the LumaRx Full-Body Hair Removal System near bathtubs, showers, basins, or other vessels containing water. Failure to follow these precautions can cause electrocution or electric shock. Contact Customer Service for assistance (page 50).
- **DO NOT** operate the LumaRx Full-Body Hair Removal System with an unapproved power cord, or if it has a damaged plug or cord. Keep the LumaRx Full-Body Hair Removal System cord away from heated surfaces (such as stoves, fireplaces, radiators or furnaces). Failure to follow these precautions can cause electrocution or electric shock. Contact Customer Service for assistance (page 50).
- **DO NOT** open the LumaRx Full-Body Hair Removal System (except to replace the removable TREATMENT CAP and/or LAMP CARTRIDGE) or attempt to repair your LumaRx Full-Body Hair Removal System as it contains no user-serviceable parts. Trying to open the LumaRx Full-Body Hair Removal System may also damage the device and will void your warranty. Please contact Customer Service (page 50) if you have a broken or damaged device. Failure to follow these precautions can cause electrocution or electric shock.
- Hair removal by lasers or intense pulsed light (IPL) can cause increased hair growth in some individuals.
- If you've had professional or clinical or self-administered laser or IPL or electrolysis hair removal treatments, wait at least 6 months before using the LumaRx Full-Body Hair Removal System. If you've waxed, plucked or epilated hair in the treatment area, wait at least 6 weeks before using the LumaRx Full-Body Hair Removal System. The device uses the hair under your skin to disable your hair at the root, so it is not effective if the hair has been pulled out.

Risks

▲ It is important to use the LumaRx Full-Body Hair Removal System properly. Make sure that you follow all the warnings (page 16), and that you do not use the device if it is not right for you. If you do not understand any of the warnings, or have any questions, please contact Customer Service (page 50).

If you use the LumaRx Full-Body Hair Removal System as directed, the risks are low:

- Your skin may become red within 24 hours after treatment. This redness generally clears up within a few hours. If the redness does not go away within 2–3 days, contact Customer Service (see page 50).
- You may develop temporary blisters, although these are unlikely for people with lighter skin tones. If you do develop a blister, stop using the LumaRx Full-Body Hair Removal System, clean the blister and cover it with a bandage and change daily to help prevent the blister from becoming infected. If you need more information, contact Customer Service (see page 50).
- In rare cases, an outbreak of cold sores may develop in the area(s) treated that have a history of cold sores. The cold sores should heal within 5 to 7 days. If the cold sores do not heal after 7 days, contact Customer Service (see page 50).
- In rare cases, some users might see an increase in the number of hairs growing. This increase is usually temporary and due to hair growth cycles. However, if after beginning use of the LumaRx Full-Body Hair Removal System you believe that the number of hairs regrowing in the treated area has increased since you started treatment, stop using the LumaRx Full-Body Hair Removal System and contact Customer Service (see page 50).
- Skin swelling, especially around the hair follicle where you may notice the skin is slightly higher than surrounding areas. If the swelling doesn't go away within

Risks

- 2–3 days, stop using the LumaRx Full-Body Hair Removal System and contact Customer Service (see page 50).
- Skin discoloration. If you notice that your Skin Tone in the treatment area has changed, stop using the LumaRx Full-Body Hair Removal System and contact Customer Service (see page 50).
 - Eye Injury. Flashing any laser or IPL device near or around the eyes could cause serious and permanent eye injury. This is very unlikely as long as the device is used as indicated.