


### Cleaning your HairMax LaserComb

- Keep your LaserComb clean.
- Turn off the LaserComb by pressing and holding the power button (lower button) until the light turns off.
- Unplug the power cord from the device and from the electrical wall outlet.
- Clean the laser window with a soft, slightly moistened cloth.
- To clean the comb teeth, gently remove any hairs or dust with a slightly moistened towel. Don't use any acetone, soap or detergents, only warm water.

- To clean the teeth more thoroughly, slide the teeth towards the top of the LaserComb and gently pull out. Clean with soap and water.
- To replace the teeth, align the guides and slide down until they lock into place. If you need to replace the teeth, you can order online at [hairmax.com](http://hairmax.com) or call Customer Service toll-free at 1-866-527-3726.

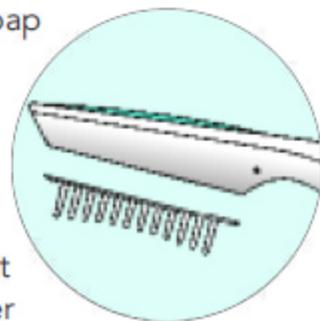

### Electrical Safety

Basic safety precautions should be observed when using this, or any other electrical device, especially when children are present.

**⚠ WARNING:** Keep away from water. Contact with water may result in an electrical shock.