

# glo<sup>®</sup>

GLO Brilliant™  
Personal Teeth Whitening Device



[www.gloscience.com](http://www.gloscience.com)

## We're GLO Science™

GLO Science brings together breakthrough science and intuitive design in beauty and oral care. Committed to bringing groundbreaking innovations in teeth whitening and oral care to consumers at home.

Invented by Dr. Jonathan B. Levine, visionary product inventor, oral health expert and NYC aesthetic dentist, your GLO Brilliant™ Device is guaranteed to get your smile white.

It's time to smile, put it on, press power and GLO!

## What's inside




## Quick Start to GLO Brilliant™ Whitening



### Get Ready to GLO®

Start with a fully-charged control and place it around your neck. Remove mouthpiece from case and firmly attach it into the side of the control.

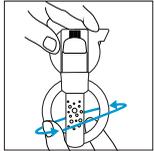

### 1 Apply Lip Treatment & Open G-Vial


Start with a fully-charged control and place it around your neck. Remove mouthpiece from case and firmly attach it into the side of the control.



### 2 Apply Whitening Gel

In a circular motion apply a thin layer of gel to the top and bottom teeth, fronts only. Gel will foam.

**Tip:** Use one G-Vial for up to four whitening applications in the session. Apply GLO Lip Treatment before each application.

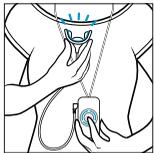

### 3 Press GLO & Whiten

Press the GLO button on the control, and your mouthpiece lights up a brilliant blue. Place mouthpiece in your mouth. The device gently warms up with each application. After eight minutes, your device will blink and automatically shuts off.

**Tip:** It is ok to swallow – the gel stays on your teeth. Repeat for a total of four back-to-back applications

## How to GLO®

You are now hands-free and have eight minutes to multitask. Apply a facial mask, answer emails, pack a lunch!

After each 8-minute application, spit out excess saliva – no need to rinse your mouth. Rinse mouthpiece in warm water and begin your next application. Repeat for a total of four back-to-back applications using one G-Vial – it's just 32 minutes to GLO and a more brilliant smile!

When complete, rinse mouthpiece and place back in its case. Whiten for five consecutive days, and depending on starting shade, you can whiten up to 10 days in a row. Your GLO Device will keep your smile white all year long. Purchase extra G-Vials to keep your GLO.

Get ready to smile more – it's never been easier!



## Charging your GLO Brilliant™ Device

It takes just 4 hours to fully charge.

- Attach small end of USB into the bottom of your control and connect the other end into the power adapter. Plug power adapter into electric outlet.
- Device blinks a brilliant blue when charging. Charge is complete when all three bands are solid and no longer blinking.
- When device is fully charged, remove GLO® Control from power cord. Lights will automatically turn off. You are now ready to GLO!



1/3 charged



2/3 charged



fully charged

**Tip:** Your GLO Control will only turn on when charged, connected to the mouthpiece, and the GLO button is pressed. Traveling and on the go? Charge your device on your computer.

Please register your device once opened. Visit [gloscience.com/registration](https://gloscience.com/registration)

## How to use

### To start Gloing and Whitening

- One whitening session consists of four, 8-minute back-to-back applications. A total of 32 minutes per day.
- Use one G-Vial per session. Each G-Vial contains enough gel for up to four, 8-minute applications.
- Use GLO® for five consecutive days. Results will vary depending on your starting shade and genetics. You can whiten up to 10 days if needed.

### For optimal results:

- Brush teeth before you use GLO Brilliant™.
- Whiten for five consecutive days. For optimal results, do not skip days in between whitening.
- Avoid food or drinks that would stain a white shirt! Use Everyday GLO pen to keep stains from setting.

### Keep your smile white and bright

- You own your GLO Device now, and can keep your smile white all year long. Touch up with one G-Vial (up to four applications) per month or as needed.

## GLO Brilliant™ Whitening Mouthpiece

The patented GLO Brilliant device uses GLO (Guided Light Optic) technology, combining elements of heat and light built into a closed system mouthpiece that accelerates the whitening process by accelerating the breakdown of hydrogen peroxide. The mouthpiece is a universal size that whitens the top and bottom teeth at the same time.

Smiles just got smart!



### Tip: Share GLO with the family

Simply purchase an additional mouthpiece for each person and you can share the same GLO Control.

## GLO Brilliant™ Whitening Gels (G-Vials)

Our unique GLO Brilliant G-Vial is designed to be used with the GLO Brilliant Device.

The G-Vial contains a highly viscous formula that stays in place when applied to your teeth. The unique brush tip targets teeth by avoiding gums and soft tissue. The G-Vial is hermetically sealed to keep the formula fresh until opened.

Each G-Vial can be used for up to four applications. Use one G-Vial for each day that you whiten, completing four consecutive, 8-minute applications. Your set of five G-Vials is good for up to 20 complete applications.



### Tip: Be prepared

Purchase an extra set of G-Vials and use to keep your smile white all year long.

## GLO® Tips & FAQs from Dr. Levine

### GENERAL WHITENING TIPS:

#### **Do I brush my teeth before I GLO?**

Yes, it is recommended that you brush your teeth before using GLO.

#### **How white can I get my smile?**

The average shade change among participants after 5 days of use was 5 shades on the upper teeth and 4.8 shades on the lower teeth.

#### **Will the brightness and speed of everyone's results be the same?**

Everyone's teeth are unique and has their own shade of white, the GLO experience will be different for everyone so take the time treat your smile individually!

#### **Will my teeth look too white and unnatural?**

No. The good news is you control the shade of white and will get the ideal shade just for you.

#### **Do I need to apply the GLO Lip Treatment before each application?**

Yes, Dr. Levine recommends it, as it will help keep your lips moist

### WHO CAN USE:

#### **Can I use GLO Brilliant with braces?**

No, but you can start as soon as they come off. In the meantime, you can try our other GLO products like Everyday GLO and Whitening Toothpaste.

#### **What causes tooth staining and discoloration?**

The most common reasons for dental discoloration could be genetics, aging, medications, colored drinks such as coffee, wine, cola and tea or illness.

#### **Can a teenager use the GLO Brilliant™ Device?**

Yes. For teenagers, age 13 to 16, complete only two 8-minute applications per day for a maximum of five days. Maintain with one 8-minute application per month. Teens 17 and older should follow adult instructions.

#### **How do I find a GLO Science™ Professional Dental Practice?**

Visit [glosciencepro.com](http://glosciencepro.com) and find a GLO Doctor near you.

### ORAL CARE TIPS:

1. Have a professional cleaning 3x a year
2. Brush 2x a day, for two minutes
3. Use a soft-bristle brush. Hold the toothbrush at 45-degree angle to the gum and brush in a circular motion
4. Floss every day
5. Brush and scrape your tongue frequently
6. Use alcohol-free mouthwash 1x a day
7. Drink water frequently to wash away leftover food from the teeth and tissues. Water also helps prevent dry mouth

## GLO® Tips & FAQs from Dr. Levine

### WHITENING GEL & G-VIAL TIPS:

**Can I reuse an opened G-Vial whitener?**  
Yes, each G-Vial contains enough whitening gel for four separate sessions. The gel remains active up to seven days with the cap closed in the refrigerator.

### Where is it best to store my extra G-Vials?

It is best to keep them in a cool dry place.

### Is it OK to swallow the gel while whitening?

Yes, the gel is okay if swallowed.

### What do I do if I experience sensitivity?

Be sure to apply GLO Lip Treatment generously to the inside and outside of gums and reapply in between applications. If sensitivity persists, visit your dental professional for a check-up as it could indicate periodontal problems.

### Can I use other whitening gels with the GLO Brilliant™ Device?

No, our G-Vials are specifically formulated and tested to work ONLY when used with GLO components. Other gels & products may cause sensitivity and unknown/ untested results.

If you are having difficulty with your device or it is not working, please call us at 855-456-5976.

### MOUTHPIECE & CONTROL TIPS:

### Will the light from the GLO Brilliant Mouthpiece or GLO Control affect my skin?

No, there is no UV light. It is white light from LEDs.

### How long will my control battery last if fully charged?

1 band of light = one application  
2 bands lit = three applications  
3 bands lit = up to nine applications

### What do I do if I need to pause my session?

First press the GLO button to turn off. Do not unplug mouthpiece while the control is lit.

### Can I whiten while charging the control?

No. For your safety, the GLO Control will not turn on while charging.

### What do I do if my mouthpiece won't light up?

Your control will ONLY work, when the mouthpiece is attached, check to be sure the cord is securely plugged into your charged GLO Control.



It's time to GLO!  
#getgloing

## GLO® One Year Limited Warranty

GLO Science™ LLC commits to the original retail purchaser that we will replace your GLO Brilliant™ Personal Teeth Whitening Device and/or components within one year from the date of purchase, for any manufacturing defect in the GLO Brilliant Personal Teeth Whitening Device or components. (Components include: the GLO Mouthpiece and Case, GLO Control, GLO Dock, Lanyard and USB/Power Adapter, GLO Brilliant Whitening Gels and GLO Lip Care).

If the product you have shipped back is a model no longer being made and is unavailable at that time, we will happily replace the product with the most recent model.

To process a return, please call 855-GLO-5976 (855-456-5976). We will provide you with a valid Return Authorization Number, which must be accompanied with your return together with a valid itemized receipt for the returned goods. To ship back your Device, for your convenience, we will provide you with a UPS prepaid postage label mailed directly to you for the return of your product. If you choose to return your Device or components and not use the UPS label we provide, we will not be responsible for the product being lost or damaged.

### Our GLO Limited Warranty does not cover any of the following

- Defects, damage or other conditions caused, in whole or in part, by misuse, abuse, neglect, negligence, alteration, accident, freight damage, improper repair, improper maintenance, improper, abnormal or excessive use and/or failure to obtain a replacement in a timely manner, as prescribed by GLO Science LLC.
- Defects or damages that are a result of an act of nature (such as a flood or earthquake) or acts of war or terrorism.
- Damage resulting from improper use of parts or accessories that are not compatible with GLO Brilliant Personal Teeth Whitening System and/or components, including but not limited to the breakage or fracture of the GLO Brilliant Personal Teeth Whitening System and/or components by any reason other than by manufacturing defect.
- Ordinary wear and tear on GLO Brilliant Personal Teeth Whitening Device and/or components.
- Defects or damages as a result of any modification of GLO Brilliant Personal Teeth Whitening Device and/or components.
- Cosmetic damage, including scratches or dents.
- Accessories or parts that are not those of GLO Science.
- Products purchased, new or used, from unauthorized distributors or through unauthorized distribution channels.

- Damage resulting from use of GLO Brilliant Personal Teeth Whitening Device with a voltage source, other than the voltage source supplied with the GLO Brilliant Personal Teeth Whitening Device.
- GLO Science's only obligation under this warranty is the replacement of defective material that comprises the GLO Brilliant Mouthpiece and the GLO Brilliant Personal Teeth Whitening Device. GLO Science shall not be liable and hereby disclaims any direct, special, indirect, incidental, exemplary or consequential damages or delays or losses under any other legal theory, including but not limited to, loss of use, loss of revenue, or loss of profits.

### Warranty Disclaimer

This warranty is GLO Science's only warranty and is in lieu of all other warranties, express or implied. GLO Science makes no implied warranties of any kind, including any implied warranties or merchantability or fitness for a particular purpose.

All of us at GLO Science are committed to making your GLO experience as positive as it can be. Simply contact us anytime at 855-GLO-5976 (855-456-5976) and we will be happy to assist you. Visit [www.gloscience.com](http://www.gloscience.com) for news and updates on GLO Science.

### GLO Brilliant Whitening Gels (G-Vials)

Net wt 5 x 2.2 ml, 5 x 0.08 oz

**Ingredients:** Hydrogen Peroxide, Distilled Water, Polyethylene Glycol 400, Glycerin, Triethanolamine, Sodium Saccharin, Versene EDTA, Carbopol, Flavor.

### Warnings:

- If irritation (such as redness, swelling, soreness) of the gums or the mouth occurs, discontinue use and consult a dentist.
- Products containing peroxides are not recommended for use by children under 12 years of age.
- Use for periods of longer than 14 days is to be only under the supervision of a dentist.
- Avoid contact of the product with the eye.
- Avoid direct contact of the active surface of the tooth whitening product with the gums and/or salivary flow.

### GLO® Lip Treatment

Net wt .30 oz, 8.50 g

**Ingredients:** Ricinus Communis (Castor) Seed Oil, Beeswax/Cire d'abeille, Mentha Piperita (peppermint) leaf, Stearalkonium Hectorite, Butyrospermum Parkii (Shea) Butter\*, Vanilla Planifolia Fruit extract, Menthol, Alcohol, Aloe Barbadosensis Leaf Juice\*, Copernicia Cerifera (Carnauba) Wax, Tocopherol, Ammonium Glycyrhizate, Ubiquinone. \*Organic.

**Cautions:** Consult with a dental professional before use if you have untreated dental cavities or any form of gum disease. If gum irritation occurs, stop using the product and contact your dental professional. Avoid getting the gel on fabric or clothing as it may lead to discoloration.

smile more

laugh more

#getgloing

[www.gloscience.com](http://www.gloscience.com)