

quick start guide


Please make sure SKINCLINICAL® reverse is fully charged before proceeding with these steps. To fully charge the device, attach the USB wall mount charger and plug it into the power outlet until “Charge Status” indicator on the bottom of the product changes from flashing BLUE color to solid BLUE color. Disconnect charger for treatment.



1. Press the on/off button to power-on.



2. The blue and green lights on the handle will illuminate and the LED lens will glow red, indicating the device is ready for use.

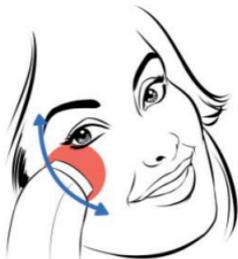

3. Place the device lens in contact with your skin and glide along treatment area until it beeps (3 minutes). Repeat for each treatment area.



4. Press the on/off button to power-off.