

# HAIRMAX®

The science of hair growth.®


## GETTING STARTED | CHARGING

### Getting Started

Read this user manual completely before you start to use your device. The instructions in this user manual are intended for individuals with a minimum of an eighth grade education.

HairMax Laser Devices are safe for home use. The patient is the intended operator of this equipment.

When you receive your device, take all the parts out of the case. Make sure your device is clean and dry before beginning the treatment.

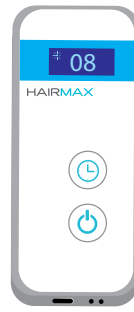

Controller/Battery Pack

Your HairMax Laser Device comes with a rechargeable Controller/Battery Pack.

You must charge the Controller/Battery Pack for 6-8 hours before initial use.

**ONLY UTILIZE THE PROVIDED CONTROLLER/BATTERY PACK WITH YOUR HAIRMAX LASER DEVICE. PLUGGING INTO ANY OTHER POWER SOURCE MAY DAMAGE THE DEVICE AND VOID YOUR WARRANTY.**

### Charging




Charge your HairMax Laser Device by plugging the small end of the Power Cord into the socket on the bottom of the Controller/Battery Pack (Figure A).

Plug the other end of the Power Cord into the USB socket of the Power Supply (Figure B). Plug the Power Supply into an electrical wall outlet.

The indicator lights on the side of the Controller/Battery Pack will flash during the charging process. Four steady blue Indicator Lights indicate your device is fully charged (Figure C).

When your HairMax Laser Device is turned on, the indicator lights will illuminate to display the battery level. As the battery charge diminishes, the Indicator Lights will turn off in a sequence to indicate the remaining charge level of the battery:

- 4 lights** = full
- 3 lights** = 75%
- 2 lights** = 50% (time to recharge)
- 1 light** = 25%

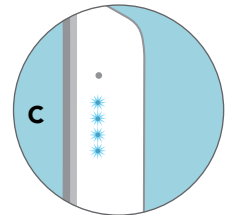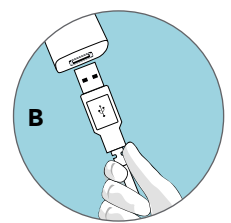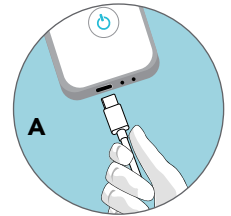

Controller/Battery Pack Side View

## PROPER USE

### Before Starting Your Treatment:

1. Make sure the battery is FULLY CHARGED per the instructions on page 11.
2. Test your HairMax Laser Device for sensitivity (see page 6 of this User Manual).

### Proper Use

Your HairMax Laser Device should be used three times per week, on non-consecutive days. For example, Monday, Wednesday, Friday.

You can use your HairMax Laser Device for 15 minutes in the default, Continuous-On Mode or 30 minutes in the Pulse Mode. See page 13 for instructions on how to select mode. Choose only one mode per treatment. Either mode provides the identical laser therapy treatment.

We recommend you use the default 15-minute, Continuous-On Mode. If you feel uncomfortable warmth in this mode, switch to the 30-minute Pulse Mode. Since the lasers are switching on and off in the Pulse Mode, less warmth is generated.

When charging is complete, disconnect the Power Cord from the Controller/Battery Pack and store it for future charging.

Connect the Controller Cable to the USB socket on the bottom of the Controller/Battery Pack, the same socket you used to charge the battery (Figure D).

Plug the other end of the Controller Cable into the Laser Connector coming out of the back of the HairMax Laser Device (Figure E).

The connectors on the Controller Cable are identical and can be used in either direction.

Push the cable into the Laser Connector Socket until it is securely joined.

You are now ready to start your treatment.

Continuous-On Mode



OR



Pulse Mode

### ATTENTION

The HairMax Laser Device must be plugged into the Controller/Battery Pack in order to operate.

**ONLY UTILIZE THE PROVIDED CONTROLLER/BATTERY PACK WITH YOUR HAIRMAX LASER DEVICE. PLUGGING INTO ANY OTHER POWER SOURCE MAY DAMAGE THE DEVICE AND VOID YOUR WARRANTY.**

### TREATMENT INSTRUCTIONS

Place the device on your head in a comfortable relaxed position (Figure F).

Press the Bottom Button (Figure G). You will hear a beep and the display will illuminate (Figure H).

By pressing the Top Button while the display is still flashing (within the first four seconds), the device will toggle between the default 15-minute Continuous-On Mode and the 30-minute Pulse Mode.

We recommend using the default 15-minute Continuous-On Mode.

The display will flash for 4 seconds and the Laser Indicator Light on the Controller/Battery Pack will illuminate (Figure H).

The Laser Indicator Light on the Controller/Battery Pack (Figure H) will indicate your lasers are on and functioning properly.

Once your treatment is completed (either 15 minutes or 30 minutes), the Controller/Battery Pack will beep twice and turn off.

If you choose to power off your device before the treatment is complete, press the Bottom Button until you hear the double beep (Figure G).

Once your treatment is completed, disconnect the Controller Cable from the Controller/Battery Pack. Then charge the Controller/Battery Pack as indicated on page 11 so it will be ready for your next treatment.

**NOTE – your HairMax Laser Device cannot be used while being charged.**

D

