


# GETTING TO KNOW YOUR TRIA AGE-DEFYING LASER

## WHAT'S INCLUDED:


- TRIA AGE-DEFYING LASER
- CHARGING CRADLE
- POWER SUPPLY AND CORD

- INSTRUCTIONS FOR USE
- QUICK START GUIDE

## IMPORTANT FEATURES:

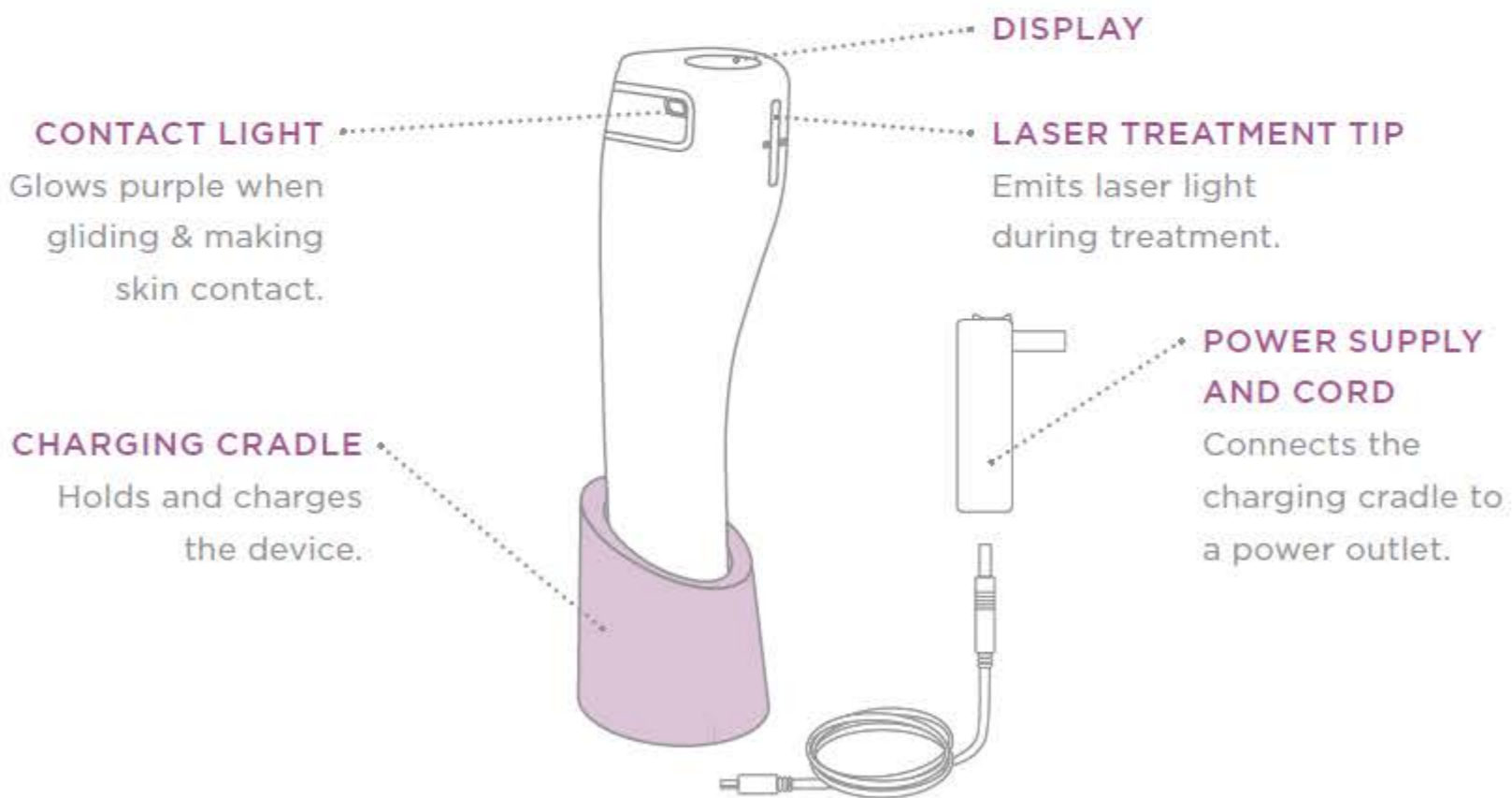

TRIA AGE-DEFYING LASER

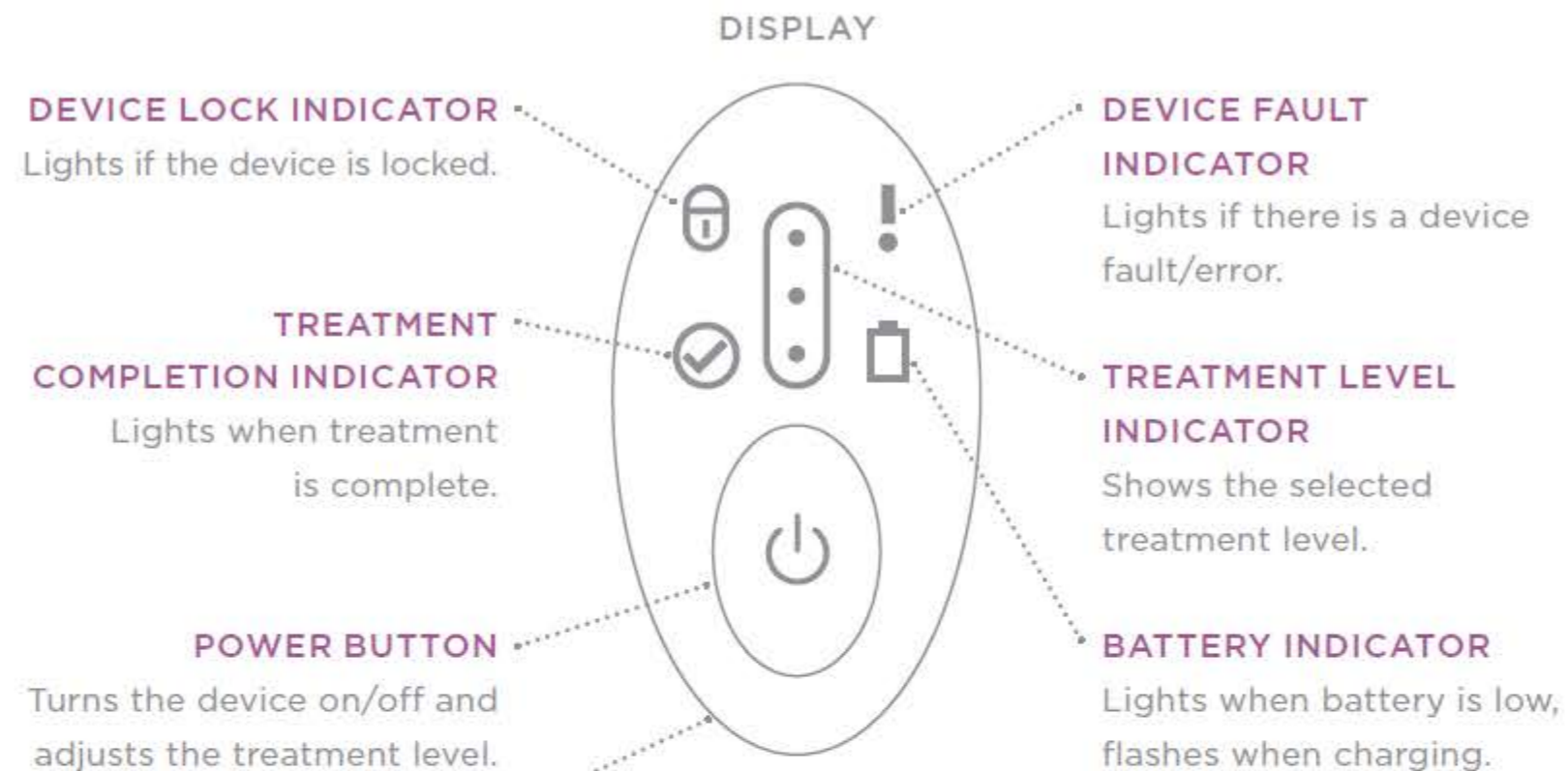

## LOCK FEATURE

When the device is locked it will not power on. The blue padlock icon will light when you press the power button.

**To lock the device:** First, turn on your Tria Age-Defying Laser. Then press the power button for several seconds until the padlock icon illuminates on the display.

**To unlock the device:** Press and hold the power button until you hear four beeps. Then release your finger from the power button.