

# TROPHY SKIN

DISCOVER YOUR BEST SKIN



## MiniMD

PORTABLE HOME MICRODERMABRASION SYSTEM

### User Guide

REV 05/2016

# MiniMD

PORTABLE HOME MICRODERMABRASION SYSTEM

We're proud to offer modern technology backed by old school customer care. We're only happy if you're happy! Have a question about your device? Contact us anytime. We'd love to hear from you! Your feedback is what helps make the MiniMD the best device it can be. Customer Care: [www.trophyskin.com/pages/contact](http://www.trophyskin.com/pages/contact)

## About MiniMD

Big benefits in a little package! The MiniMD is a compact and convenient tool designed to deliver salon-quality microdermabrasion treatments. This lightweight and portable glow-getter uses a high tech combo of diamond exfoliation and gentle suction to buff away dead skin cells, revealing the radiant, glowing skin that sits just below the surface. MiniMD effectively polishes and revitalizes skin for an instantly softer, smoother complexion. It's small enough to fit in your purse, travel case, or gym bag so you can treat your complexion to microdermabrasion anywhere, anytime!

## With regular use, MiniMD:

- Deeply exfoliates
- Invigorates the complexion
  - Smooths and softens
- Polishes and revitalizes skin
- Buffs away dead skin cells
- Improves the look of skin radiance

The Art of Microdermabrasion: [www.trophyskin.com/video](http://www.trophyskin.com/video)  
Register your warranty at: [www.trophyskin.com/warranty](http://www.trophyskin.com/warranty)

## Table of Contents

4	Getting to Know MiniMD
5	Precautions & Care
6	Usage Guide
8	After Treatment
9	Troubleshooting
10	Warranty
11	Optional Accessories

## Getting to Know MiniMD



## Precautions & Care

### PLEASE READ ALL INSTRUCTIONS BEFORE USE

When using electrical products, follow basic safety precautions including the following:

#### **IMPORTANT: SKIN SENSITIVITY TEST**

Before using MiniMD for the 1st time, we advise you perform a skin sensitivity test on the inside of your forearm. Do this by using High level and glide the diamond tip across the entire surface of your inner forearm using the tips in instructions on Page 6. After 15 minutes, if you notice any hairline scratches or uncomfortable irritation, please contact us for a fine tip for sensitive skin.

1. Always unplug this product immediately after use
2. Keep away from all liquids
3. Do not use while bathing and do not place product where it can fall into a tub or sink
4. Do not reach for a product that has fallen into water. Unplug immediately.
5. Do not disassemble the RejuvadermMD for any reason. Contact us if you need assistance.
6. Keep a black wool filter in place at all times to prevent foreign substances from entering the motor.
7. Handle with care. Dropping the device may damage the frame or motor.
8. Do not use the diamond tip around the eye area, areas that contain inflammatory acne, or areas that have rashes or wounds.
9. Stop using the device immediately if you develop any abnormal reactions such as itching, uncomfortable inflammation, or scratches.
10. Do not use the RejuvadermMD device for commercial purposes.

## Complete Usage Guide

### PREPARE FOR USE



#### STEP 1:

Thoroughly wash your face with a foaming cleanser of your choice and pat dry.



#### STEP 2:

Connect the power adaptor to a power outlet, then plug the adaptor into the bottom of the MiniMD.



#### STEP 3:

Press the Power button to turn on the device.



#### STEP 4:

Press the Power button again to change from suction level 1 to suction level 2. We recommend starting with level 1 for the first few uses until skin becomes accustomed to treatment.

## Treat Your Skin

#### STEP 1:


Brace the skin with your free hand, and place the Diamond Tip flat against the skin.

#### STEP 2:

Using short, even strokes, gently move the Diamond Tip across the skin, ensuring that the tip remains flush with skin.

#### STEP 3:

Treat areas such as the face, neck, chest, arms, hands, or feet 2-3 times per week using this method. You can increase the suction level over time if needed.

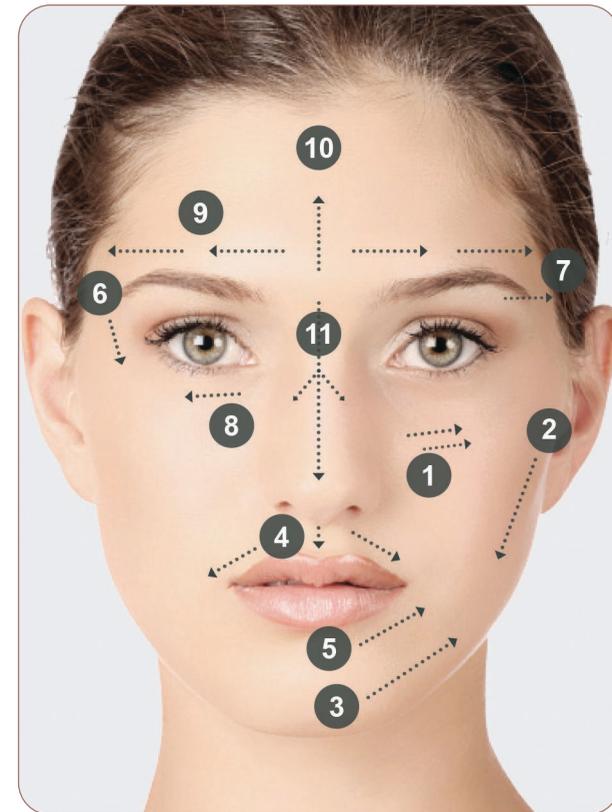

#### CHEEKS:

Stroke with the Diamond Tip following Arrows 1-3, while bracing skin with other hand.

#### AROUND LIPS:

Follow Arrows 4 and 5. Do not touch your lips with the Diamond Tip.

#### AROUND EYES:

Stroke downwards between eye and ear (Arrow 6), outwards on temple (Arrow 7), and below eye (Arrow 8). Do not stroke eyelid and stay outside of the eye orbit.

#### FOREHEAD:

Stroke outwards, and work your way upwards with horizontal strokes (Arrows 9 & 10).

#### NOSE:

Cover entire nose as per Arrow 11.

## After Treatment

### STEP 1:

Gently slide off the Diamond Tip, and discard the black wool filter.

### STEP 2:

Replace with a new filter and push the Diamond Tip back into place.

### STEP 3:

Wash your face with cold water to remove remaining debris and gently pat dry.



### Step 4:

Apply our Argan Stem Cell Serum (not included) or another high quality anti-aging serum immediately after each treatment.

You may experience some redness and sensitivity after your treatments, which will usually subside within a few hours. Use 2-3 times per week to help keep your skin looking smooth and radiant.

### PRO TIPS FROM OUR TEAM

- Don't stroke more than 2-3 times in the same spot.
  - Apply the diamond tip in slow, smooth strokes to ensure that the tip remains flush with the skin.
- Gently brush the Diamond Tip with an unused toothbrush after each use to remove dead skin cells and debris.
  - Apply SPF after use as skin is extra sensitive to sun exposure post-treatment.

## Troubleshooting Guide

### My device won't turn on

Make sure the power adaptor is plugged in securely at the back of the unit, and plugged into a working power outlet. Press and hold the Power button for a few seconds. If the light does not come on, please contact customer service.



### My device has lost suction

Check the Diamond Tip closely for any hairline cracks, punctures, or spaces. Remove the Diamond Tip and ensure that both rubber O-rings are in place. If you notice cracks or missing O-rings, please contact us for a replacement.



Hairline Crack

### My Diamond Tip isn't coarse enough

The Diamond Tip will accumulate skin cell build-up over time, making it feel smoother. Buff the Diamond Tip with an unused toothbrush regularly to keep it clear of skin debris. The Diamond Tip should be replaced every 3-6 months. A replacement can be purchased on our website.



## Satisfaction Guarantee

Trophy Skin is committed to 100% customer satisfaction. If for any reason you are unsatisfied with the results of your MiniMD device, you may initiate the return process by sending an email to [sales@trophyskin.com](mailto:sales@trophyskin.com) within 30 days of receiving your device. The following terms will apply to be eligible for our full, "no questions asked" refund (less S&H) policy:

- MiniMD is in new condition.
- You have followed the treatment instructions outlined in this manual.
- Device must be shipped back at your cost with tracking number, using the shipping method of your choice.
- All included accessories must be returned and must be in new condition.

## Warranty

Our first priority is your happiness! Your MiniMD is warranted to be free from defect in material and workmanship for a period of 6 months under normal use after its original purchase date. This warranty extends only to the original retail purchaser and only when purchased through the Trophy Skin website or approved 3rd party affiliated sites.

If your MiniMD should become defective during the warranty period, contact us via email at [sales@trophyskin.com](mailto:sales@trophyskin.com) or call **866-546-4664** for repair or replacement. We reserve the right to replace a defective product with the most comparable product currently available. We will work with you to ensure that you are satisfied.

To view our full Limited Warranty (and to register your warranty) please visit: [www.trophyskin.com/warranty](http://www.trophyskin.com/warranty) or call our office at **866-546-4664** during normal business hours.

## Disclaimer

Trophy Skin Is Not Engaged in the Practice of Medicine

Trophy Skin expressly disclaims any and all responsibility for any consequences, directly or indirectly, related to any use, misuse or interpretation of the website or any other information provided or omitted in connection therewith and you agree to indemnify, defend and hold Trophy Skin harmless therefrom. Trophy Skin shall not be deemed to be engaged in the corporate practice of medicine as a consequence of selling the MiniMD acne light therapy device. In addition, you specifically agree to with the following statements:

- I understand the MiniMD is a cosmetic treatment and that no medical claims are expressed or implied.
- I understand that to achieve maximum and continued results the protocol recommended by Trophy Skin should be followed.
- I understand there are no guarantees implied as to the results of MiniMD treatment, due to many variables, such as: age, skin type, skin condition, sun damage, smoking, alcohol, environmental exposures, etc. I understand that I may or may not actually see demonstrable visible results, that each case is individual.
- I understand that if I am pregnant I must consult my physician before using MiniMD or making any changes in my normal routine.

## Medical Treatment or Diagnosis Disclaimer:

Nothing contained in [www.trophyskin.com](http://www.trophyskin.com) is intended to be instructional for medical diagnosis or treatment. The information presented by [www.trophyskin.com](http://www.trophyskin.com) and/or Trophy Skin should NOT be considered complete, nor should it or take the place of a consultation with a physician or competent healthcare professional for medical diagnosis and/or treatment, but is provided as information only. The MiniMD listed and sold on this site should ONLY be used after consultation with a dermatologist, general physician or other qualified health care provider.

Information obtained on this site related to acne treatment is NOT meant to be exhaustive and obviously cannot cover all situations and conditions. Again, we urge you to seek consultation or advice of your physician or other qualified health care provider promptly with any health care related questions. You should never disregard medical advice or delay in seeking it because of something you have read on this web site or others.

The information contained in [www.trophyskin.com](http://www.trophyskin.com) is compiled from a variety of information providers. The web site [www.trophyskin.com](http://www.trophyskin.com) does not directly or indirectly practice medicine or dispense medical services. Better Business Bureau (BBB) Binding Arbitration and Dispute Resolution. Any complications or disputes that arise out of the purchase or use of the MiniMD device will be resolved through BBB Binding Arbitration program. As a member of the BBB, Trophy Skin elects that all disputes with customers will be resolved under this program.

## Optional Accessories



### Body Diamond Tip

Our diamond body tip makes quick work for full body treatments. It has a coarse, wider surface area to treat larger body parts.



### Precision Diamond Tip

The precision tip has a fine grain and small surface that is perfect for those hard to reach, delicate and/or angular areas of the face.



### Fine Diamond Tip

The fine coarseness diamond tip is suitable for those with sensitive skin.



### Argan Stem Cell Serum

Custom-formulated to be used after microdermabrasion treatments.



### Replacement Filters

Use only 1 filter per treatment and purchase a 100 Pack of Replacement Filters as needed.



### Hydrating Facial Mist

Hydrating Facial Mist provides instant hydration for your face and neck after microdermabrasion and light therapy treatments.

# TROPHY SKIN

DISCOVER YOUR BEST SKIN

**TROPHY SKIN CUSTOMER SERVICE**  
469.698.6270 **TOLL FREE** 866.546.4664

[SALES@TROPHYSKIN.COM](mailto:SALES@TROPHYSKIN.COM)

TROPHY SKIN  
4372 KENMARE TRAIL  
FRISCO, TX 75034

Distributed by:  
Ideal Sourcing Ltd, PE1 5WG, UK  
Made in China

