
# LEGXERCISE

*Professional*

Contador de Pasos & Control Remoto

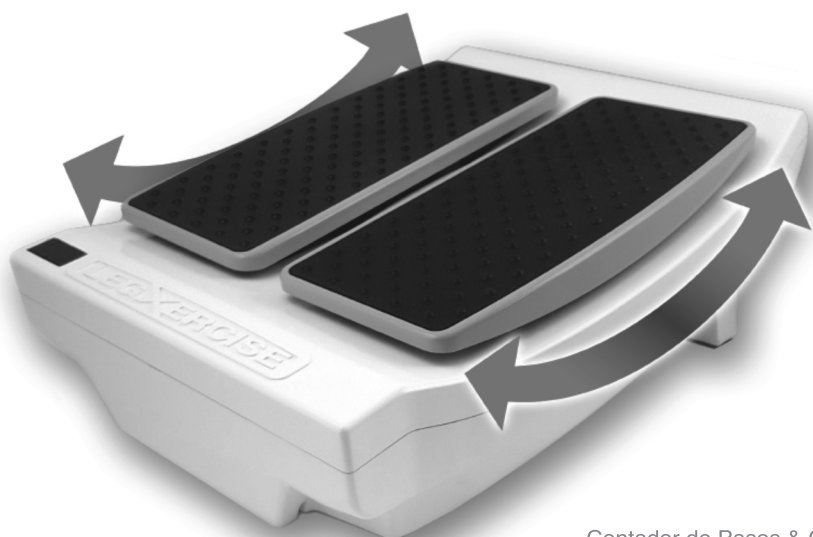
## Instructions :: Instrucciones

2 AAA 1.5V  
Batteries Required

To open

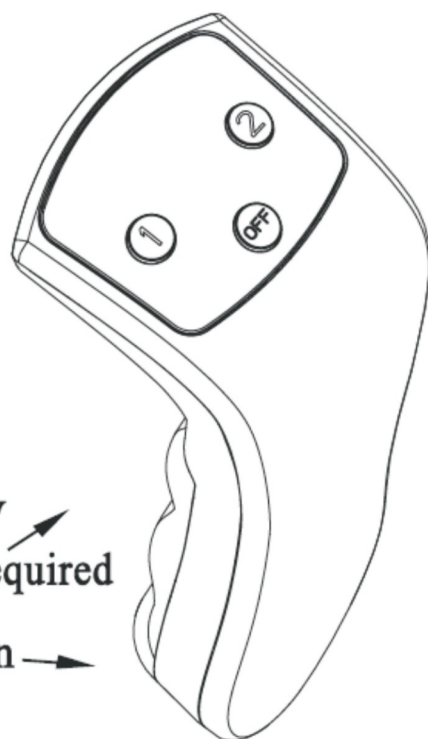
Legxercise is a Revolutionary Exercise Machine that will keep your legs in movement while sitting, provides Passive Exercise that can be performed while sitting, at work, watching TV, or reading. Beneficial Continuous Motorized Movement, without any physical effort on your part.

*Walking Simulator Motion Tech*, maximum overall leg movement, power for 70 steps per minute, new comfort remote control, increased knee movement, plus effective ankle flex, new digital step counter, 4,200 steps (2.66 km) (1.65 miles) per hour.

## The complete system to help maintain your physical health

The back and forth movement promotes rest and relaxation for your feet.


Ideal for all those people that spend a lot of time sitting or standing.

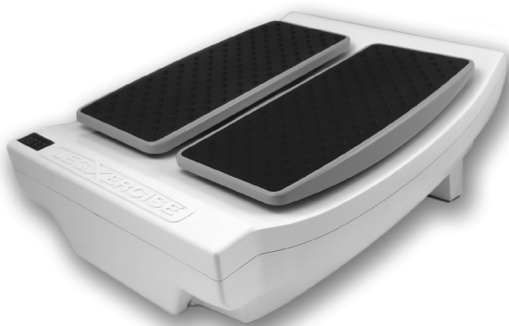An easy-to-use product from the comfort of your home or your workplace.

Helps to improve blood circulation.

### Select the Desired Speed

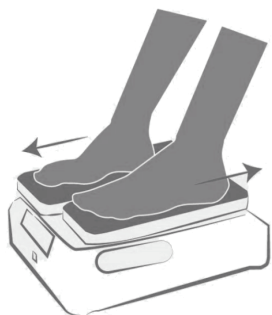
**(1) Normal Speed:** Recommended for continuous and normal use.


**(2) Fast Speed:** Recommended when you want to take advantage of all of its power for faster results.


## Instructions for Use

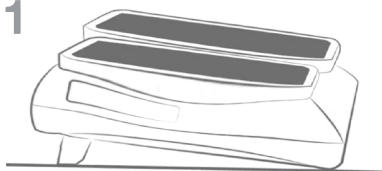**1**


Unit's location: must be flat on the floor.  
It has two rear legs. Pull out if necessary.

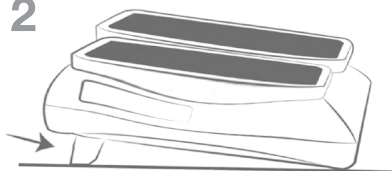**2**

We recommend using some level of inclination, it helps with posture.

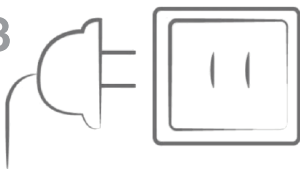**3**

Plug the adapter into the nearest electrical outlet. Plug the other end of the adapter into the front of the machine.

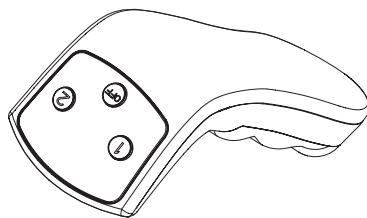**4**

Select the desired operating speed:  
Normal Speed 1 and Fast Speed 2.

**5**

Place your feet on top of both footrest,  
**DO NOT STAND UP.** There is no need to  
apply pressure or put forth any effort, just  
sit back and relax, allow Legxercise to do  
all the work for you. It will work without  
any physical effort on your part.

**6**

Use Legxercise Professional as many  
hours a day as you feel comfortable while  
seated, long use will show better results  
due to a healthy leg circulation.

**WARNINGS:**

- This unit should not be exposed to water or splashes of water. Objects that contain liquids like a flower vase, coffee cup, etc., should not be near the unit.
- If you use a power strip, keep the strip away from humidity or damaged cables.
- Keep fingers away from moving parts at all times.
- **KEEP OUT OF REACH OF CHILDREN.**
- To clean it, turn the unit off and unplug it. Then wipe clean with a dry cloth.
- Unplug the unit while not in use.

Legxercise was designed to be used while seated only. **NEVER STAND UP ON THE MACHINE**, as the weight of your body will damage the machine's internal mechanism and its functions.

**How to get rid of electrical and electronic old appliances:**

1. If the symbol with a crossed container is displayed, this means it complies with Directive 2002/96/CE.
2. All electrical or electronic appliances should be disposed of properly at designated local places only.
3. The correct disposal and treatment of damaged goods contributes to avoiding potential risks to the public and the environment.
4. To obtain more information on how to properly dispose of your electrical and/or electronic devices, contact your local authorities or the retail place of purchase.

**DANGER**

KEEP HANDS & FEET AWAY  
FROM IN BETWEEN PEDALS  
& MACHINE AT ALL TIMES.  
KEEP OUT OF REACH OF  
CHILDREN.

**⚠ CAUTION**

Do NOT  
stand on machine.

**Height Adjustable Levers**

Pull out levers under  
machine to adjust  
your position.

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