


# Quick Setup Guide



## Features



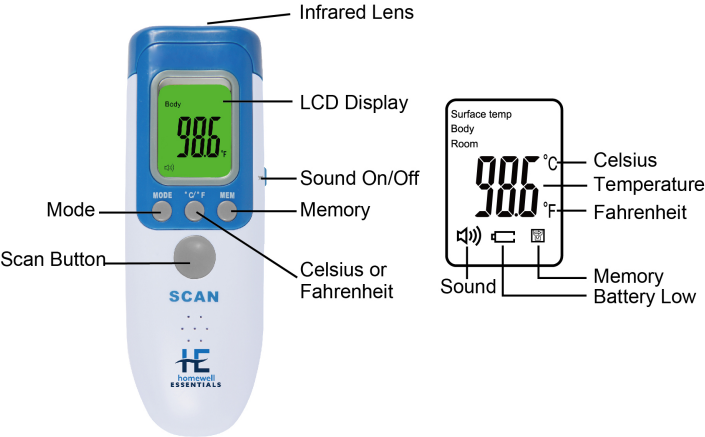

## Start Here

1. Install the included batteries. Make sure the terminals (+/-) are correct.
2. Wait for 10 to 15 minutes for the thermometer to acclimate to the ambient temperature before use.
3. Press the Scan Button, point towards the forehead from a distance of 1.2 to 2 inches away and press Scan Button again to take the temperature.
4. Understanding the numbers.



## Troubleshooting

If you are experiencing problems while using your thermometer, please refer to the chart below to help resolve the problems. If the problem still exists, please contact our customer service.

Thermometer doesn't turn on.	Check the battery installation. Please make sure the terminals are aligned correctly.
The body temperature seems inaccurate.	It could be caused by the dramatical change of ambient temperature. Please wait for 15 to 20 minutes before using it.  Or it could be caused by hair or perspiration or air draft. Please move the hair, wipe the perspiration or move to the location where there is no air blowing.
The temperature seems very off.	Please check if it is set in the right temperature mode. The body temperature can only be measured in Body mode. The surface temperature mode measures the external area temperature.
Showing Lo on the screen. 	It does not mean the battery is low. It could be caused by pointing the thermometer to the wrong place when it is the Body mode.  In Body mode, the display shows Lo when the temperature is lower than 89.6°F(32°C)
Showing Hi on the screen. 	It could be caused when by pointing the thermometer to the heat source when it is the Body mode.  In Body mode, the display shows Hi when the temperature is higher than 109.2°F(42.9°C)

## 5 Year Limited Warranty

Consult the user guide for the full warranty policy.